




Installation and Assembly: Medical Antimicrobial Base Shelf

Model: ACC315-AB, ACC315-AW

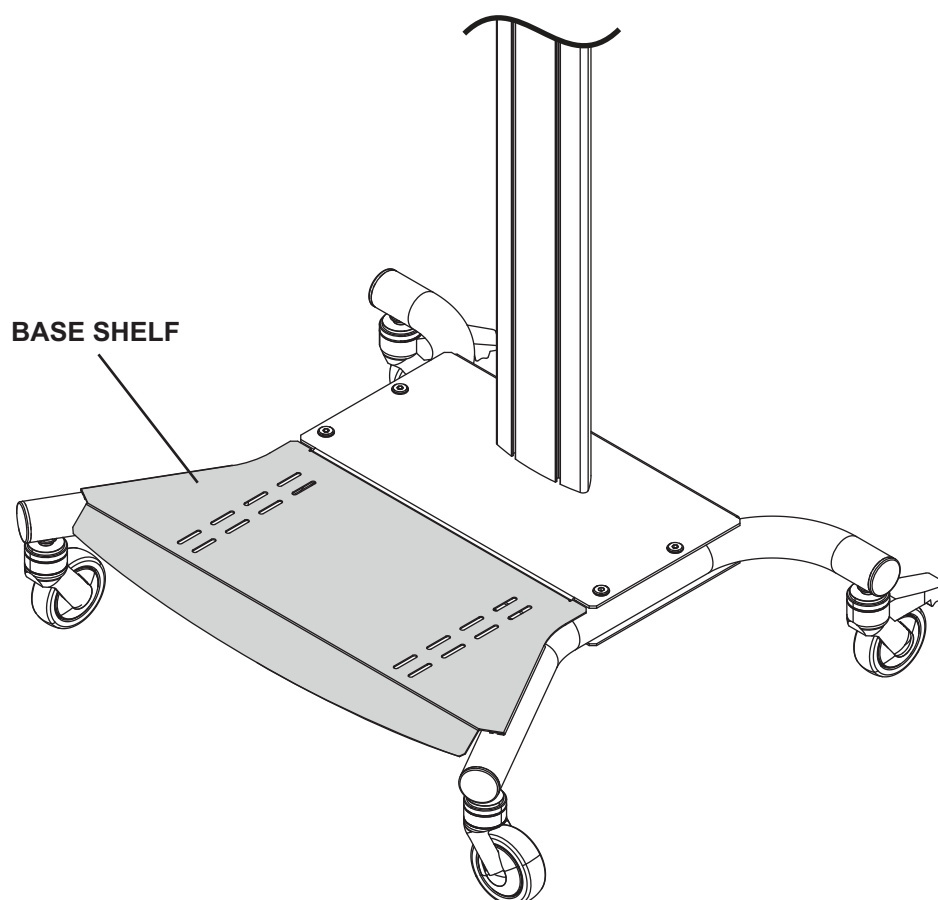



Features:

- Treated with Agion® antimicrobial finish in black or white designed to assist with infection control.


Before you begin, make sure all parts shown are included with your product.

Parts List

Description	Qty.	ACC315-AB	ACC315-AW
		Part #	Part #
A extended base shelf	1	009-S1226	009-S2226
B M5 x 12 mm phillips pan head screw	2	520-1027	520-1027
C washer	2	540-9400	540-9400
D split lock washer	2	540-1035	540-1035

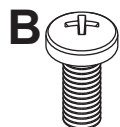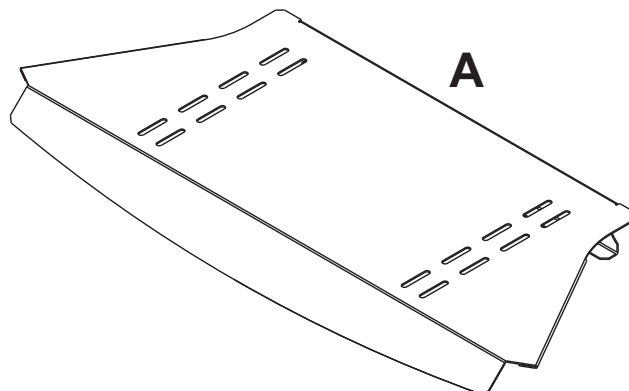

- 1** Attach extended base shelf (A) to cart base as shown, using two M5 x 12 mm phillips screws (B), two washers (C) and two split lock washers (D).

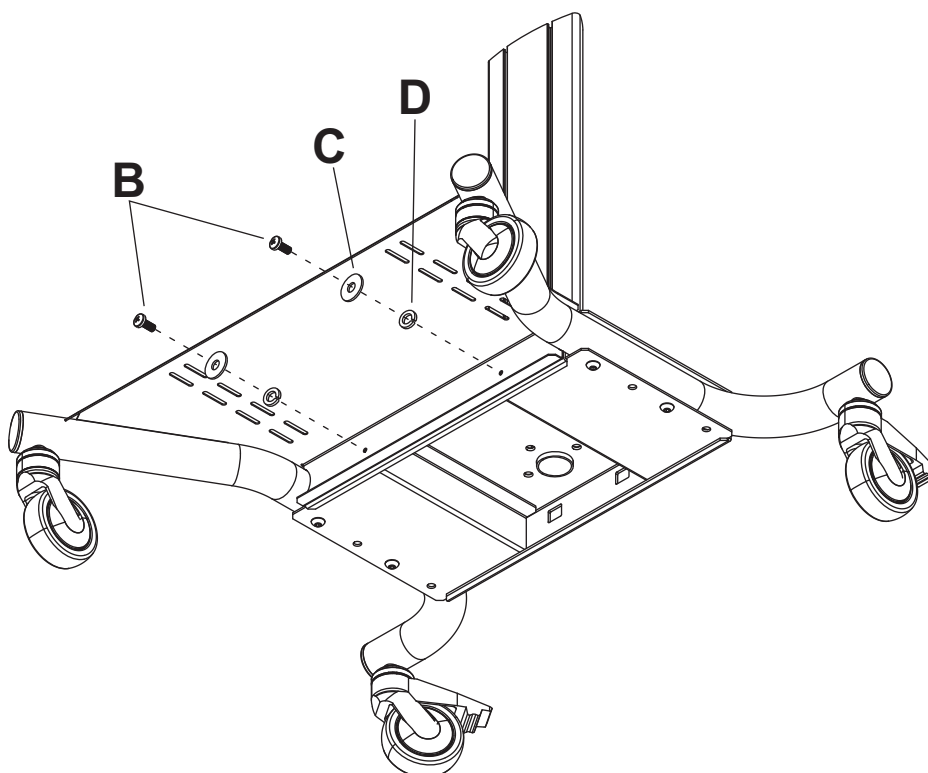

Recommended Cleaning Guidelines

The following procedures are not guaranteed to control infection. An infection control administrator or epidemiologist should be consulted regarding cleaning procedures and processes.

Most painted components will withstand cleaning by commonly used, diluted, non-abrasive solutions such as quaternary ammonia compounds, ammonia enzyme cleaners, and bleach or alcohol solutions. However, it is recommended that any surface be tested in an inconspicuous area prior to determine the propensity for discoloration.

- Pen or permanent and dry erase markers can be removed with 91% isopropyl alcohol and a soft cloth.
- Iodine stains can be removed with commonly used cleaners and a soft cloth.
- Never use steel wool or other abrasive materials that will damage the surface finish.
- Do not use Acetone, mineral spirits, abrasive cleansers, paint thinner or any other harsh or toxic chemicals to clean your product. These solvents will damage the surface finish.

WARNING

- To avoid risk of electric shock, do not expose electrical components to water, cleaning solutions or other potentially corrosive liquids or substances.
- Do not use flammable cleaners on product surfaces due to close proximity of electrical power and equipment.